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## Recommended Usage Method

STEP 1: After cleansing the face, apply an adequate amount of gel or essence. Press and hold the "LED/☺ button" to turn the device On/Off, press the "MODE button" to select the mode, press the "LEVEL button" to select the intensity, and press the "LED/☺ button" to select the color light. There are a total of 5 modes, 4 levels of intensity (excluding the eye mode), and 3 types of color lights. Please choose according to your needs.

STEP 2: When the button is pressed correctly, a "Dee" sound will be emitted. Press the "Left Side/Right Side/Forehead/Eyes button" and follow the recommended care sequence to switch between them, and then slide to use. When one side is finished, a "Dee" sound will be heard as a prompt to remind you to switch to another side.

STEP 3: When the device emits a long "Dee~" sound, it indicates that the current mode of treatment has ended. Clean the remaining gel from your face and clean the treatment head.

TIPS: During use, move the device in close contact with the skin and replenish gel or essence as needed to achieve the best results and ensure the normal operation of the device. When the indicator light comes on, please charge the device promptly.


For detailed instructions on how to use the product, you can scan the QR code to watch the instructional video.

### **Note:**

Do not use on the eyelids, skin around the eye socket, eyeball, or throat. Be sure to use water-soluble skincare products. Do not use oily or exfoliating skincare products, as this may cause device malfunctions or skin problems.

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
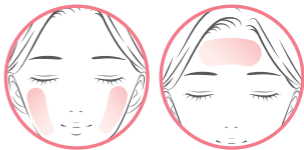
RF / EIP / EMS / COMBO mode：  
Please use the device by sliding  
in the following order.

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EYE mode: When using the EYE mode,  
please follow the recommended order  
of use and take caution not to apply it to  
the shaded areas.

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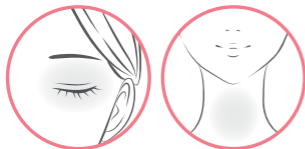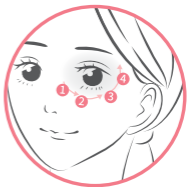Please avoid using the device  
on the eyelids, skin around the  
eye socket, eyeball, or throat.

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## Skincare Advice

Mode	Recommended Skincare Products*	Time	Intensity	Frequency in Use(Ref)**
<b>RF</b> Wrinkles Reducton & Skin Firming	Specialized Gel / Essence	Est. 10 mins	4 Levels	2-3 times/week
<b>EIP</b> Effective Penetration	Specialized Gel / Essence	Est. 10 mins	4 Levels	2-3 times/week
<b>EMS</b> Contour Lifting	Specialized Gel / Essence	Est. 10 mins	4 Levels	everyday
<b>COMBO</b> Quick Comprehensive Care	Specialized Gel / Essence	Est. 5 mins	1 Level	2-3 times/week
<b>EYE</b> Eye Care	Specialized Gel / Eye Essence	Est. 5 mins	1 Level	2-3 times/week

Please use the color LED function as needed in combination with the above modes.  
Red LED: Activates and revitalizes the skin, improves blood circulation.  
Orange LED: Brightening and reducing dark spots, enhances skin immunity.  
Blue LED: Anti-inflammatory and oil control, improves acne and pimple problems.

\* Please be sure to use water-soluble skincare products. Do not use oily or exfoliating skincare products as they may cause mechanical failures or skin problems.

\*\* This information is for reference only and can be adjusted according to individual needs.