# Timeless PA-BD01

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User Manual

Timeless

Paige Luminous™

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1. Preparation for Use	46
2. Product Details	57
3. Functions	58
4. Operation	59
5. Skincare Advice	64
6. Treatment Sequence	65
7. Treatment Aftercare	66
8. Product Specification	68
9. Warranty	69
10. Disclaimer	72

# Table of contents

Thank you for choosing our products! To ensure proper and safe use, please take the time to read this manual carefully and follow the instructions provided. After reading, please keep the manual in a safe place for future reference. If you plan on transferring the device to another user, be sure to provide them with the manual as well, so they can also familiarize themselves with its usage.

## Preparation for Use

### Safety Precautions:

- Warning: Not following the instructions could result in negative consequences
  - Please read the "Safety Precautions" carefully before using the device
  - Proper usage requires following the steps correctly.
  - If you notice any abnormalities, stop using the device immediately
  - Using the device in an improper condition can result in harm to the body
  - To prevent personal injury or damage to property, the below diagrams included in this manual clearly indicate the potential extent of injury and loss
-  **Danger:** An accident that could result in death or serious injury.
-  **Warning:** An accident that may cause death or serious injury.

 **Caution:** An accident that could result in minor personal injury or property damage, including damage to household furniture, livestock, or pets.

 **Prohibited:** Indicates prohibited matters.

### **Danger**

**Do not use in situations that may cause medical electronic equipment to malfunction.**

If you are using any of the medical electronic equipment listed below, please do not use this device. Otherwise, incorrect operation of the medical electronic equipment may cause accidents and physical discomfort."

- Cardiac pacemaker and other implanted electronic medical instruments
- life support equipment such as artificial heart and lung machines
- wearable electronic medical instruments such as electrocardiography monitors and infusion sets

### **Warning**

**This device should not be used by the following individuals as it may result in accidents or bodily injury:**

- Those who are currently receiving medical treatment or have fever symptoms as it may interfere with their medical condition
- Visceral disease patients

- Tuberculosis patients
- Patients with heart disease, abnormal blood pressure, and malignant tumors
- Patients with hemophilia, bleeding disorders, diabetes, facial neuralgia
- Individuals with weakened immunity due to diseases such as AIDS or who are using immunosuppressants
- Individuals undergoing dental treatment or have had dental fillings
- Individuals with abnormal perception, sensation of pain or inability to sense temperature
- Individuals with keloids, facial melanosis, photosensitivity, or dyspigmentation
- Individuals suffering from telangiectasias caused by prolonged use of steroid hormone drugs or liver dysfunction
- Individuals with allergies, hereditary dermatitis, atopic dermatitis, skin diseases, sensitive or fragile skin
- Individuals with a history of epilepsy or who experience discomfort from exposure to flickering lights
- Individuals with sequelae of the head or eyes due to traffic accidents or other similar incidents
- Individuals who have had face or eye surgery, including LASIK, should consult a doctor before using the device
- Individuals who cannot express themselves, including young children
- Individuals who are still developing

- Individuals with high myopia (above-6D)
- Pregnant women and women who are menstruating or lactating
- Individuals who are drunk or have a hangover
- Individuals who are taking sleeping pills
- Individuals who are extremely tired or unwell
- Individuals who have skin problems such as pigmentation and contact dermatitis
- Individuals who have had or currently have dermatitis caused by cosmetics
- Individuals who have spots or any other skin condition caused by physical stimulation
- Individuals who have skin inflammation or redness caused by sunburn
- Individuals with skin purulent, inflammation or acne
- Individuals who have recently undergone chemical peels at beauty salons
- Individuals who have experienced skin problems after using beauty equipment
- Warning: Using this beauty device repeatedly or for an extended period on the same area of skin could potentially cause damage
- Warning: Do not bypass or disable the inherent safety features of the beauty device

**Do Not Use The Device On The Following Body Parts And Areas:**

- Mucous membranes around the heart, eyeballs, eyelids, lips, and inside the mouth

- Areas of melanosis, lesion, mole or viral wart
- Areas that cause skin irritation and inflammation due to excessive sun exposure
- Areas with wounds, eczema, psoriasis, herpes, purulent acne and pimples
- Itching, fever, physical stimulation and other parts that cause pathological pigmentation
- Areas that have undergone surgery, genitals, anus, around Adam's apple, ears
- Areas that have been recently depilated or shaved
  - ✘ It is important to consult your doctor before using the device if you are receiving hospital treatment or have any medical concerns. If you are unsure about any symptoms you may have, please seek medical advice before using the device.

### **Prohibition**

- Do not use “Timeless” simultaneously with other beauty devices
- Do not use the device excessively or continuously in areas where the treatment head cannot fully fit, as it may result in burns or low-temperature burns
- Do not use the device for a prolonged period of time on thin skin or over bones
- Be sure to use the device by sliding it and avoid keeping it in one place for an extended period of time
- Avoid using on the same area continuously as it may lead to skin damage and other issues like infection, inflammation, bleeding, pigmentation, or worsen existing symptoms

- Do not use the device if it is damaged or not functioning properly
- Do not use the device for any purpose other than its intended cosmetic use
- This beauty device is intended for cosmetic use only. Do not use it for any treatment or purpose other than those specified in the user manual
- Do not look directly at the LED light
- Do not use this device on pets
- To prevent overheating, fire, or malfunction, do not expose the device to high temperatures or heat sources, and do not attempt to charge the device in heated environments or by throwing it into a fire
- Do not charge, use, or store the device in the bathroom or in areas with high humidity, or in environments with high levels of oil or smoke
- Do not place the device in areas where it may come into contact with water, such as near a sink or washstand
- Do not use or expose the device to water, as it may result in electric shock, short circuit, fire, or malfunction
- Do not use this device near a bathtub, shower, sink, or other container with water 
- Avoid touching the device with wet hands: as this may cause electric shock, short circuit, fire and malfunction
- Do not use the device near flammable items, such as alcohol, thinner, volatile oil, lighter fluid, nail polish remover, or sprays, as it may cause heat generation, fire, or malfunction

- Do not use the device in heavily dusty areas or under direct sunlight
- Do not use the device in areas with strong electromagnetic waves or magnetic fields
- Do not place in a microwave oven or high pressure container
- Avoid using in areas with sudden temperature changes, such as near air conditioners or heaters
- Immediately halt use if you detect an unusual odor or if the device becomes excessively hot. Continuing use may result in damage or fire
- Do not place heavy objects on top of the main unit, and be careful not to drop or bump it while using it to avoid any malfunctions or injuries
- Do not keep the device in places that are prone to instability, tilting, vibrations, or collisions. Avoid dropping or bumping the main unit as it can cause mechanical failure or injury
- It is not recommended for children 14 years old or younger to use this device
- Children must not use, clean, or maintain this device without supervision
- This device may not be used by individuals (including children) who have reduced physical, sensory or mental capabilities, lack experience and knowledge, unless they are under supervision or have been instructed
- This device must only be used under SELV (Safety Extra Low Voltage) conditions
- Use a power adapter with an output power no greater than 15W
- This device must not be used on areas of broken skin

- Use only the disposable and consumable materials recommended by the manufacturer
- Maximum output voltage: 32V AC.

### **Caution**

- Avoid applying excessive pressure or focusing on a single area while using the device as this can result in skin injury
- Avoid colliding the treatment head with other metal objects as it may result in heat generation, fire, failure, or damage
- Before use, inspect the treatment head for cracks, damages, deformations or other issues. Do not use the device if any of these problems are found, as it may result in skin damage or malfunction
- Do not use the device in extreme environmental conditions, such as extremely hot or cold temperatures. (The recommended ambient temperature range for using the device is between 15°C to 30°C.)
- Do not place the device on the floor or furniture when it is turned on. Turn off the power immediately after use to prevent damage to furniture or potential injury
- Avoid contact with other electronic products such as electronic blankets while the device is in use and the power is turned on. This may cause malfunction or injury
- Do not connect this device to other metallic or electronic products such as mobile phones or television sets

- Do not use this device in combination with exfoliating skin care products or masks, as this may lead to skin issues

### **Warning (About Charging)**

- Ensure that the charging port is securely inserted into the charging base, and the USB charging cable is fully plugged into the base port, when charging the device. Failure to do so may result in electric shock, short circuit, or fire
- Do not damage, bend, pull, twist, or place heavy objects on the USB charging cable, or clamp it tightly. Such actions may cause fire, electric shock, or malfunction
- When moving or picking up the device, do not pull the USB charging cable. This may cause electric shock, short circuit, fire, or poor contact
- If the USB charging cable is damaged or loose when plugged into the socket, stop using it immediately
- When connecting or disconnecting the USB charging cable, it is important to hold the front of the USB and avoid touching the metal parts
- Avoid inserting or removing the USB charging cable when your hands are wet
- Keep the USB charging cable and **base** away from water to prevent any potential hazards such as electric shock or damage to the device
- Clean the USB charging cable and **base** before charging to prevent any residue from affecting the charging process
- Unplug the USB charging cable after the device is fully charged

- To prevent heat generation or fire, avoid overloading the sockets with too many devices or appliances
- This device is intended for household use only
- If the following situations occur, please stop using the device , turn off the power and unplug the USB charging cable from the socket:
  - Thunderstorms, earthquakes, etc.
  - Wiring smoke, abnormal heat, strange smell, deformation
- If an abnormal situation occurs, immediately turn off the power and stop using the device

### **Prohibition of Disassembly**

- Do not attempt to disassemble, modify, or repair the device on your own. If repairs are needed, please consult your dealer

### **Caution on Use**

- The device may become warm during use, which is normal and poses no safety issues
- It is important to clean the treatment head before each use to prevent skin irritation or other issues
- It is important to remove any precious metal jewelry, glasses, contact lenses, or other items before using the device to prevent potential damage or irritation
  - ※ Precious metal: ear ring, necklace, ring & watch, etc.
- It is recommended to start with the lowest setting (level 1) when using the instrument for the first time or until you are comfortable with it. This will help prevent any discomfort or injury

- When using the treatment head, do not apply too much pressure on the skin, and avoid staying in one place for more than 3 seconds. To avoid skin problems, do not use on the same area repeatedly within 24 hours
- When using the device on the skin around the eyes, be sure to close your eyes to avoid discomfort or injury
- Please keep the skin moisturized between the treatment head and your skin by applying skin care products regularly
- If you experience a stinging sensation on the skin, reduce the intensity, move to a different area, and apply moisturizer before continuing use
- If you experience pain or any abnormal sensations on your skin or body while using the device, stop using it immediately and seek medical advice
- Excessive use of the device may have adverse effects on the health of your skin
- It is important to use the device in accordance with the condition of your skin
- It is recommended to apply water-soluble skincare products only after properly removing any makeup, cleanser, and debris from the face. The condition of your skin may change based on the overall health of your body
- It is important to regularly assess the condition of your skin while using the product. If you experience any symptoms such as rashes, redness, itching or other abnormalities, stop using the product and seek the advice of a specialist


## Product Details

Total Electrode Area:  
531.32mm<sup>2</sup>


24K Gold Treatment Head

The surface of the device is coated with 24K gold, which has a strong affinity with the skin. However, due to the soft nature of the metal, the surface may become scratched easily. Please note that small scratches on the surface will not affect the quality of the device and it is safe to use as normal.

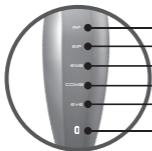Mode Indicator

Low Battery Indicator

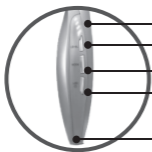Intensity Indicator /

Charging Indicator

Intensity Button

Mode Button

Charging Button / LED Button

Charging Port

## Accessories

- Charging Port & USB Type-C Charging Cable
- User Manual & Quick Start Guide
- Warranty Card & Storage Pouch
- 24K Gold Aqua Gel

## Functions

RF (Radio Frequency )

Warming the skin from the bottom using the device can help to improve the skin's firmness, elasticity, and plumpness. This can be particularly effective in addressing concerns such as wrinkles, fine lines, and sagging skin.

EIP (Electroporation Interval Pulse)

Introduces the ingredients of beauty essence into the deep layers of the skin quickly and efficiently.

EMS (Electrical Muscle Simulation)

The micro-current stimulates continuous muscle contractions, enhances the metabolism of muscle cells, and results in a lifted facial appearance.

Red LED

Promotes blood circulation, reduces inflammation, and speeds up the healing process of wounds, etc.

### Orange LED

Helps to prevent dry skin, increase skin elasticity, smooth wrinkles and stretch marks, delay the signs of aging, and other effects. It can also enhance skin immunity, reduce skin sensitivity, decompose pigments, brightening and reducing dark spots.

### Blue LED

Helps inhibit oil secretion, reduce inflammation and sterilize the skin, eliminate acne-causing bacteria, treat acne and scars, decrease the production of sebum, and improve acne and pimple problems.

### LED Fully Open


The use of red, orange, and blue LED lights is suitable for individuals with multiple skin concerns.

## Operation


When using the device for the first time, please remove the protective film (transparent film) on the treatment head.

Low battery Indicator 	
Off	Fully Charged
Red	Low Battery, need to be recharged
Flashing Red Light	Very Low Battery, recharge immediately


## Charging Using A Household Socket:

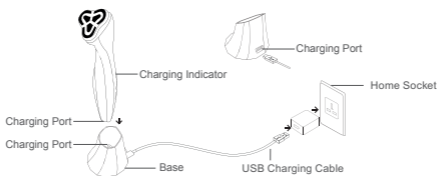
## Charge The Beauty Device Via The Charging Port:
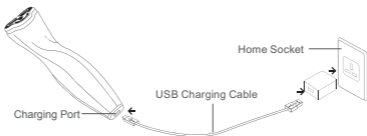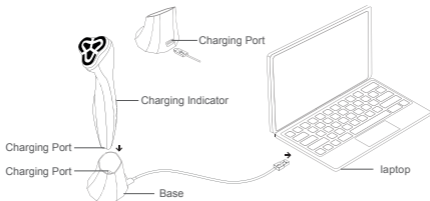
## Charging Using A Computer:
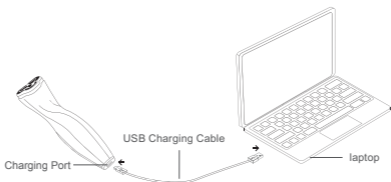
## Charging Using A Computer:

1. Please ensure that the battery is fully charged before using the device for the first time. During charging, the low battery indicator will flash red, while the charging indicator will flow. Once charging is complete,  the low battery indicator will turn off , and the charging indicator will remain lit.
2. To turn on the device, press and hold the "LED/⏻ button". When turned on for the first time, the treatment head will light up in red.
3. To select the color of the light, briefly press the "LED/⏻ button". The colors cycle through in the following order:  
Red LED → Orange LED → Blue LED →  
LED Full Brightness → off  
Caution: When the light is on, do not look directly at the source.

4. Mode selection: press briefing the "MODE button" to select the mode and cycle through the modes in the following order:  
RF → EIP → EMS → COMBO\* → EYE  
\* COMBO is a 5-minute quick comprehensive care mode that combines the benefits of RF, EIP, and EMS for a complete skin care experience.
5. Strength level selection: Press the "LEVEL button" in RF / EIP / EMS / COMBO mode to select the intensity. There are a total of 4 intensity levels to choose from.
6. The default working time for the RF / EIP / EMS modes is 10 minutes. During the first 5 minutes, treat one side of the face and after 5 minutes, there will be an audible signal indicating it is time to treat the other side. After 10 minutes, the device will automatically enter standby mode and the indicator light for the next mode will flash.
7. The preset working time for the COMBO/EYE mode is 5 minutes. Start by treating one side of the face/ eye for the first 2 minutes and 30 seconds, then switch to the other side after the beep. The device will automatically stop and enter the EYE standby state after 5 minutes of using the COMBO mode. After completing the EYE mode, the device will automatically shut down.

- By long pressing the "LED/⏻ button", the device will be shut down. The RF/EIP/EMS/COMBO mode has a memory function, and it will automatically resume at the previous intensity of that mode after restarting.

**Caution:**

- When using for the first time, it is recommended to start at level 1
- Although this device is equipped with a device that prevents low-temperature burns, for safety purposes, please be sure to use it by sliding in RF/COMBO/EYE mode
- The EYE mode is specially designed for eye care and has only one intensity level. When using this mode, it is important to keep your eyes closed

**Description of Mechanical Sound:**

- The "Dee" sound is an indication that the button has been successfully pressed.
- The "Dee Dee" sound indicates that an incorrect button press has been made.
- The "Dee Dee Dee" warning sound will be emitted in RF mode if the device stays in one position for more than 3 seconds.
- In RF, EIP, and EMS modes, a "Dee" sound will be triggered to remind switching to another side for care when half of the suggested care time has elapsed.
- At the end of the recommended care time for each mode, a long "Dee" sound will be emitted as a prompt.

## Skincare Advice

Mode	Recommended Skincare Products*	Time	Intensity	Frequency in Use (Ref)**
<b>RF</b> Wrinkle Reduction & Skin Firming	Specialized Gel / Essence	About 10 Mins	4 Levels	2-3 Times / Week
<b>EIP</b> Effective Penetration	Specialized Gel / Essence	About 10 Mins	4 Levels	2-3 Times / Week
<b>EMS</b> Contour Lifting	Specialized Gel / Essence	About 10 Mins	4 Levels	Daily
<b>COMBO</b> Quick Comprehensive Care	Specialized Gel / Essence	About 5 Mins	4 Levels	2-3 Times / Week
<b>EYE</b> Eye Care	Specialized Gel / Eye Essence	About 5 Mins	1 Level	2-3 Times / Week
<p>Please use the color light function as needed in combination with the above modes.</p> <p>Red LED: Activates and revitalizes the skin, improves blood circulation</p> <p>Orange LED: Brightening and reducing dark spots, enhances skin immunity</p> <p>Blue LED: Anti-inflammatory and oil control, improves acne and pimple problems</p>				

\* It is important to use water-soluble skin care products. Avoid using oily or exfoliating skin care products, as this may result in machine failure or skin irritation.

\*\* The recommended guidelines are for reference only and can be adjusted to meet individual needs.

## Treatment Sequences

### Sliding Motion:

The recommended sliding speed is 2cm / sec. If it feels too hot, feel free to increase the sliding speed.

### Caution:

- Ensure that the treatment head is fully in contact with the skin while using.
- Use by sliding over the skin.
- Continuously using the device on the same position may cause low-temperature burns on the skin.
- Avoid applying the treatment head on the eyelid and orbital bone.

RF / EIP / EMS /  
COMBO mode :  
Please use the device  
by sliding in the  
following order.

**EYE mode:** When using the EYE mode, please follow the recommended order of use and take caution not to apply it to the shaded areas.

Quickly learn how to use this beauty device by scanning the QR code and watching the video.

## **Treatment Aftercare**

After every use, clean the device promptly to maintain its cleanliness. Cleaning steps are as follows:

- Use a soft cloth that has been dampened with water, then wring it out, to clean the device. Avoid using hard cloths as they may scratch the surface and affect its performance
  - ※ Avoid using solvents such as alkali, alcohol, thinner, benzene, volatile oil, gasoline, or nail polish remover, abrasive powder, hot water, or medicines to clean the device.

- ※ If the device comes into contact with essence or gel, please store it after cleaning.
- After cleaning, place the main unit on the base to prevent the treatment head from being damaged by other objects

### **Avoid Storing The Beauty Device In The Following Environments:**

- Places exposed to direct sunlight, humid or dusty • Areas that are constantly exposed to direct sunlight, humidity or dust
- Places with excessively high or low temperatures, such as inside a car under the scorching sun
- Flammable areas, areas exposed to wind and rain.
- Areas with exposure to chemicals or gas and places subject to strong vibrations
- Close to heating appliances and their surroundings, as well as areas with strong magnetic fields

If the instrument has not been used for over six months, please inspect it for any abnormalities before use. (Check the LED lighting status and look for cracks, etc). Using the device in a faulty state may result in an accident.

## Product Specifications

Product name	Paige Luminous Timeless Beauty Device
Model	PA-BD01
Rated Input	5V  2A
Rated Output	About 10W
Charging Time	About 2 hrs
Usage Time	About 40 mins
Dimensions	Approx. 41mm (L) x 59mm (W) x 183mm (H)
Product Weight	Approx. 130g

\* The charging time is related to the rated input. The charging time mentioned above is based on the rated input of 5V 2A. If the rated input is lower than 5V 2A, the charging time will become longer.

## Safety Standard Certification

You can learn about the certifications obtained by this beauty device by scanning the following QR code.

## Warranty

### After Sales Service:

This beauty device comes with an original warranty of two years.


### Rechargeable Battery:

- This device contains a rechargeable battery that must be replaced by a professional
- Rechargeable batteries are considered consumables and are not covered under the warranty
- When the usage time after charging is extremely short, it may indicate that the battery life has reached its end
- In case you need to replace the rechargeable battery, please consult your dealer for assistance

### Regarding Disposal and Environmental Protection:

- When disposing of this device, please follow local laws and regulations.
- The substance contained in the built-in lithium ion rechargeable battery may be harmful to the environment, so the battery must be removed before disposing of this device.
- This device must be powered off before removing the battery, and the battery must be disposed of safely.
- Do not use other batteries to replace the battery inside this device, especially non-rechargeable batteries, as this may cause fires or electric shocks.
- Do not dispose of the battery with household waste. Ensure that the battery is disposed of at an official recycling center or collection point.

## Battery Removal Steps Illustrated:
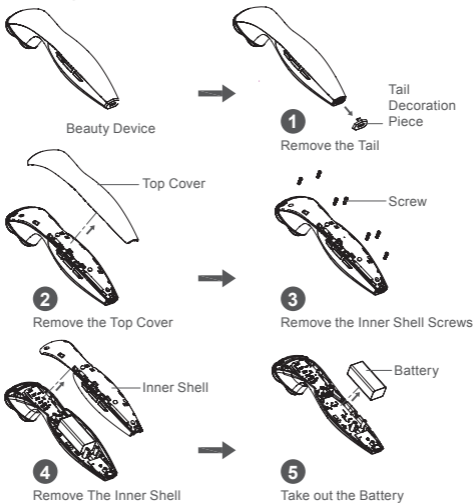
### **Contact Information:**

Hong Kong Company: Regal Silver Limited

Tel Number: 852-21866089

Address: Unit 1209, Vanta Industrial Centre, 21-33  
Tai Lin Pai Road, Kwai Chung, N.T.,  
Hong Kong

Shenzhen Company: Shenzhen Tian Hao Regal Silver  
Technology Limited

Tel number: 86-13378665715

Address: 701, Zhiyan Innovation Building, Southeast  
side of the intersection of Genyu Road and  
Tianliao Road, Tianliao Community, Yutang  
Street, Guangming District, Shenzhen,  
Guangdong, China

E-mail: [cs@regal-silver.com](mailto:cs@regal-silver.com)

To ensure a smooth process, please contact customer service by email to confirm before sending the beauty device to the address listed above.

## Disclaimer

Regal Silver Limited is committed to ensuring that all information contained in this user manual is current and accurate, but shall not be held responsible for any errors or omissions. All information is subject to change without notice.

Regal Silver Limited shall not be liable for any loss incurred due to alleged misuse or use of the product in any manner other than the usage methods outlined in this user manual, on the product packaging and leaflet. When using the Timeless beauty device, customers are responsible for understanding and complying with all product safety information.

This user manual may not be copied or translated into any language without the written authorization of Regal Silver Limited.

Regal Silver Limited makes no guarantees as to the effectiveness of Timeless beauty device. The effectiveness of Timeless beauty device will depend on the user's individual skin condition. We strongly recommend that all users consult a qualified healthcare professional before using any product(s).

The aforementioned statements have not been evaluated by the Food and Drug Administration. The products developed by Regal Silver Limited are not intended to diagnose, treat, cure, or prevent any disease. Before using any product, if you have any questions, please consult your professional healthcare advisor.

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